
Ultra-Low Temperature Freezer

User Manual

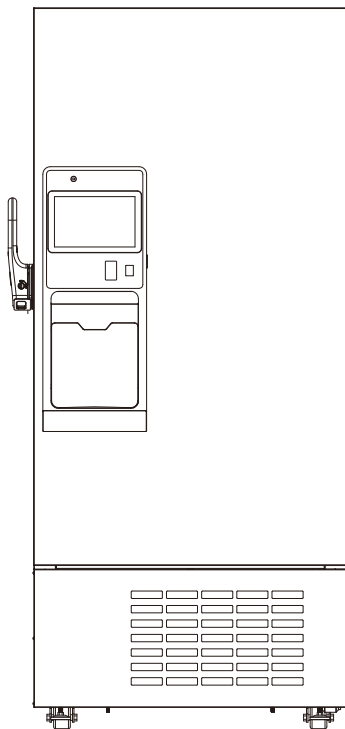MODELS:

AULT530DC | AULT780DC

BEFORE USE, PLEASE READ AND FOLLOW ALL SAFETY RULES AND OPERATING INSTRUCTIONS.

Models:	Serial Number:
----------------	-----------------------

Felix Storch, Inc.
An ISO 9001:2015 registered company
770 Garrison Ave
Bronx, New York 10474
www.summitappliance.com

Table of Contents

1.	Application Notes	1
2.	Safety Instructions	2
3.	Precautions While in Use	4
4.	Product Installation	4
	4.1 Handling and Transportation	4
	4.2 Installation Environment	4
	4.3 Installation Site	5
	4.4 Preparation Before Use	5
	4.5 First Power-on	5
	4.6 Operation after Power Failure	6
	4.7 Air Switch	6
5.	Location of Parts	7
6.	Operating Instructions	9
	6.1 Introduction of the Homepage	9
	6.2 Alarm Management	11
	6.3 Access Management	12
	6.4 Data View	14
	6.5 Data Curve	17
	6.6 Settings	17
	6.7 Message Board	21
	6.8 Login Options	22
	6.9 System Features	24
	6.10 Remote Alarm Terminal and RS-485 Interface	24
7.	Maintenance	25
	7.1 Defrost and Cleaning	25
	7.2 Equipment Discontinuation	26
	7.3 Battery Maintenance	27
8.	Troubleshooting	30
9.	Limited Warranty	32

Cutting-edge control technology:

- High-precision microcomputer temperature control system, platinum resistance temperature sensor, with temperature that can be accurately set in the set range
- LED/LCD touch screen with clear display, good man-machine interaction interface, simple operation and more accurate temperature control

Refrigeration system:

- High-efficiency compressor and EBM fan with strong power, energy conservation and high efficiency
- Double-layer foaming thermal insulation door, surrounding design of multi-patent external door system, which can effectively prevent cooling loss
- The freezer body is made of high-performance vacuum insulation material, which greatly improves the insulation effect

Humanized design:

- New power-assisted door handle, easy to open
- A special file box to make the recording simpler and more convenient
- Adjustable shelf structure, convenient for storing items
- Universal wheels and fixed feet to move and fix more conveniently

Security system:

- Keyboard lock and password protection function to prevent arbitrary adjustment of operating parameters
- Remote alarm and various alarms system(sensor failure, high temperature, low temperature, door opening, power failure, etc.) to ensure the safe storage of items in the cabinet
- Protection functions such as startup delay and shutdown interval to ensure reliable operation of products
- Door lock to prevent door opening without authorization

Main Technical Parameters

Model	Ambient temperature	Climate type	Refrigerant	Rated voltage (V~)	Rated frequency (Hz)	Temp Range (°C)	Volume (L)	Net Weight (kg)	Exterior dimensions (DxWxH) (mm)
AULT530DC	16~32°C	N	Mixture Gas	110-230	50/60	-40~-86°C	530	276	1041 x 930 x 1947 41" x 36 5/8" x 76 5/8"
AULT780DC	16~32°C	N	Mixture Gas	110-230	50/60	-40~-86°C	780	331	1025 x 1205 x 1955 40 3/8" x 47 1/2" x 77"

* The foaming material of this product is cyclopentane.

Accessories

Operation Manual	Certificate of Conformity	Key (s)	Defrosting shovel
1	1	2	1

** The included accessories may vary; please refer to the actual items delivered.

1. Application Notes

Thank you for choosing our product. For safe, convenient operation and proper maintenance of the equipment, please read this Operation Manual carefully before use and retain it for future reference.

Operators may copy selected sections of this manual for internal use only, such as for training or emergency procedures. These sections are clearly identified in the manual's table of contents.

The manufacturer assumes no responsibility for equipment damage or malfunction resulting from failure to operate the equipment in accordance with the instructions and methods specified in this manual.

There are currently no known undesirable side effects associated with the use of the ultra-low temperature freezer when operated as instructed.

Due to continuous product improvements, the functions described in this manual may differ from those of the product you have purchased. In such cases, refer to the actual product and its installed features.

- Please carefully read the Attention and Safety Precautions in Section 2: Safety Instructions before operating the equipment.
- During transportation and use, avoid violent vibration or impact. Protect the freezer from rain. Store the unit in a clean, well-ventilated environment with humidity not exceeding 80% and free from corrosive gases.
- Maintenance and repair of the equipment must be performed only by the manufacturer or an agent authorized by the manufacturer.
- If the operator encounters a situation not addressed in this manual, contact the manufacturer or an authorized service agent for proper guidance.
- Failure to operate the equipment in accordance with the specifications and instructions in this manual may result in damage to the equipment.
- Use only accessories provided or approved by the manufacturer. The manufacturer assumes no responsibility for adverse consequences resulting from the use of unapproved accessories.
- Inspect and maintain the equipment regularly to ensure proper and reliable operation.
- The Ultra-Low Temperature Freezer is intended for cryopreservation of plasma, cells, tissues, and similar materials. It is suitable for use in blood centers, hospitals, epidemic prevention centers, research institutes, and similar facilities.
- Do not use mechanical devices or other means to accelerate the defrosting process, except those recommended by the manufacturer.
- Any serious incident involving the device must be reported to the manufacturer and to the applicable regulatory authority in the country of use.

Warning:

After each use, be sure to dry the water stains around the sealing strip to prevent freezing. Wear protective gloves when storing items in the equipment, and beware of frostbite!

Tips

- Always use protective equipment correctly (including clothes, gloves, goggles, etc.)
- Always practice good hygiene habits
- Each personnel is obliged to be responsible for his or her own safety

2. Safety Instructions

When using this product for the first time, please review the meanings of the following warning symbols and carefully read all safety precautions to ensure safe and proper operation of the equipment.

 Warning! Failure to observe the precautions may result in serious personal injury or death.	 Attention! Failure to observe the precautions may result in personal injury or equipment failure and related property losses
 Protective conductor terminal Must be grounded	 Risk of Explosion This sign indicates the risk of explosion when using volatile and explosive chemicals.
 Warning: Pinch point	 Warning: Low temperature
 Warning: Fire hazard	

 Warning: Failure to observe the precautions may result in serious personal injury or death.
 If a CO ₂ backup system is adopted, the equipment use site must be provided with a good ventilation environment, and the high concentration of CO ₂ in the air will endanger the life safety of personnel. In the case of poor ventilation, please consider other measures to ensure a safe working environment.
 Do not touch the equipment with wet hands to avoid electric shock. Equipment shall not be grounded through gas pipeline, water supply pipeline, telephone line, or lightning rod, which may easily cause electric shock.
 Be sure to install the equipment on a solid and flat ground and take due care to prevent tipping over. If the ground is not solid enough or the installation location is not appropriate, it may cause the equipment to fall over and cause equipment damage or personal injury.
 Please handle the power cord carefully to avoid short circuits or open circuits. Please turn off the power before pulling out the power plug. Hold the power plug carefully and pull it out. Do not pull the wires of the power plug. Otherwise, it may cause electric shock or fire due to the short circuit.
 Please insert the power plug into the outlet tightly to ensure firmness and reliability, so as to avoid electric leakage. After installation, the power plug must be within reach, so that the power plug can be unplugged in time in case of emergency.
 Separate special outlets must be used and grounded reliably. The cross-sectional area of the copper conductor in the wall connected with the outlet must be more than 4mm ² . Do not lengthen the power cord without authorization to avoid heating or fire.
 Do not use the power supply that is not specified in the equipment design, so as to avoid overheating, short circuit and other faults. For example, connecting 110V rated voltage products to 220V power supply may cause faults such as overheating and equipment burning. For detailed input voltage AC please refer to Specification (Rated voltage $\pm 10\%$). If the operating voltage is too low or too high, a suitable automatic voltage stabilizer must be installed for cooperative
 Place the equipment on stable surface to avoid shaking.
 Do not place the equipment in a dangerous area, and do not operate the equipment near flammable items to prevent explosion or fire accidents.

	Do not place the equipment in areas exposed to the sun or rain, so as to prevent danger such as short circuit or overheating.
	Do not tilt or lay the equipment sideways, and do not impact the equipment body. Refrigeration systems are installed in the equipment, which is easy to be damaged by tilt or impact.
	Please place the equipment in a dry and dust-free environment to avoid risks such as overheating, and short circuit.
	In case of unexpected sound, smell, smoke, etc. when the power is turned on, please unplug the power in time and contact the manufacturer or supplier.
	Please place the equipment in a dry and ventilated environment, and ensure that the equipment vents and instrument surfaces are not blocked or shielded by walls or other objects. Do not use it in a poorly ventilated environment to prevent damage caused by heat released by equipment.
	It is forbidden to disassemble and modify this equipment without authorization, so as to avoid potential safety hazards. In this case, we will not bear any responsibility for quality accidents.
	It is forbidden to put inflammable and explosive dangerous goods, strong corrosive acids, alkalis and other items unsuitable for the equipment in the equipment.
	Metal objects such as nails or iron wires shall not be inserted into any aperture and gap or any outlet of the equipment, otherwise electric shock or injury may be caused due to accidental contact between the above objects and moving parts.
	In order to ensure the normal operation and ventilation and heat dissipation of the equipment, the back, left and right sides of the freezer body should be at least 30cm away from the wall, and the air inlet and air outlet must not.
	This equipment must be connected to a ground wire, and the equipment conductor terminal must be protected.

	Note: Failure to observe the precautions may result in personal injury or equipment failure and related property losses.
	It is forbidden to store living animals, flowers or other items with strict temperature requirements in the equipment.
	The temperature in the freezer is very low during operation. Do not touch the liner surface covering plate and frozen objects with your hands without wearing protective gear or frostbite may occur.
	Hold the handle and close the door to avoid pinching your fingers. When the equipment is not used for a long time, unplug it and pack it for storage.
	When restarting the equipment after power failure or power off, check the equipment settings first, otherwise the stored items may be damaged due to the change of settings.
	The equipment can be used for item preservation, not as production equipment!
	When handling the equipment, please be careful not to tip over the equipment, so as to prevent equipment damage or personal injury.
	When handling, it shall be lifted from the bottom, with the inclined plane not be greater than 45°, and it shall be handled with care. Use the equipment in safe areas. Improper use may cause harm to human health or environment.
	Pay attention to checking the filter screen in time and clean it if necessary, so as to avoid the temperature rise or failure of the equipment caused by poor ventilation due to dust jam.

3. Precautions While in Use

- When the equipment is running, the contact part between the front part of the equipment and the outer door may generate heat. This is not a fault: to prevent condensation around the casing, a heating anti-condensation pipe is installed in the equipment.
- Before placing items into the equipment, ensure that the temperature in the freezer has reached the set value first, and then put the articles in batches. Do not load more than one-third (1/3) of the freezer's internal volume at a time, as excessive loading may cause a temporary temperature rise.
- The freezer cabinet is equipped with an access port for routing a test line during testing. After the test line has been routed, the access port must be resealed with thermal insulation material. Failure to reseat the access port may prevent the freezer from reaching the set temperature and may cause condensation to form around the exterior of the access port.
- The temperature displayed on the control panel represents the temperature measured by the sensor inside the freezer. When the unit is first started, there may be a difference between the displayed temperature and the actual temperature at the center of the chamber. As the freezer reaches stable operation, the displayed temperature will gradually align with the actual internal temperature.
- Clean the equipment using a diluted neutral cleaner only. Do not use brushes, acidic cleaners, gasoline, soap powder, polishing agents, or hot water, as these may damage the painted surfaces and plastic or rubber components. Do not wipe plastic or rubber parts with volatile solvents such as gasoline.
- After a period of operation, frost may accumulate on the inner walls and inner door of the freezer. Excessive frost buildup can reduce insulation performance and increase power consumption. When the frost layer reaches approximately 5 mm, use the supplied defrosting tool to remove the frost.
- Before defrosting, remove all stored items and place them in a suitable storage environment to prevent damage caused by temperature rise during the defrosting process.
- Cooling coils are located behind and along the sides of the inner walls. Do not use sharp objects such as knives, ice picks, or screwdrivers to remove frost. Avoid scratching the inner walls during defrosting, as this may result in equipment damage or malfunction.
- If the equipment will not be used for an extended period, disconnect the power supply.

4. Product Installation

4.1 Handling and Transportation

- The product is heavy, the handling process requires the use of forklifts or hydraulic trucks inserted underneath the bottom bracket, so as not to cause damage to personnel and products.
- Avoid touching and scratching the preservation box during handling. When handling, the maximum tilt angle cannot exceed 45° (the angle between the horizontal direction), so as not to cause refrigeration system failure.
- Packaging removal: the product is packaged in a wooden bottom bracket, column packaging, and cardboard box.

4.2 Installation Environment

- Ambient temperature: 16°C ~ 32°C, the most ideal temperature is 18°C ~ 25°C, and the air conditioning system shall be used when necessary.
- Relative humidity: $\leq 80\%RH$.
- Ensure there is no strong vibration and corrosive gas around.
- Avoid areas with large amounts of dust.
- Avoid rocking or shaking the equipment.
- Elevation of the working position of the equipment: less than 2000m.
- Indoor use, pollution degree 2, and over voltage category II.
- For detailed input voltage AC refer to Specification (Rated voltage $\pm 10\%$).
- Avoid areas with direct sunlight, other cold and heat sources and strong electromagnetic interference, which will affect the normal operation of the control system and directly damage the system in severe cases.

Attention: The ambient temperature has a significant impact on equipment performance. If the specified environmental requirements are not met, the equipment may not operate properly. Ensure the environment meets the required conditions before using the equipment. The equipment is intended for intermittent operation.

4.3 Installation Site

To ensure proper operation and optimal performance, the installation site must meet the following requirements:

- Do not install the equipment in a narrow or enclosed space. The room doorway must be wide and tall enough to allow normal movement of the equipment. Adequate access must be maintained to prevent maintenance difficulties in the event of failure, which could otherwise lead to damage to stored items due to delayed repairs.
- The installation surface must be solid, level, non-combustible, and capable of supporting the full operating weight of the equipment.
- Install the equipment in a well-ventilated area and protect it from direct sunlight.
- Each unit must be connected to a dedicated power outlet. Ensure the power plug is securely connected to the outlet.
- Verify the rated operating voltage before use. In locations with unstable power supply, use an appropriate voltage stabilizer rated for the equipment load to ensure the input voltage meets installation requirements.
- The equipment must be properly grounded. If the power outlet includes a grounding connection, verify that grounding is effective before use. If the outlet is not grounded, grounding must be installed by qualified engineers or technicians before operating the equipment.

4.4 Preparation Before Use

1. Remove the outer packaging of all products (packaging requirements: The product is wrapped with PE wrap film, the internal buffer material is honeycomb cardboard and foamed polystyrene, and the outer packaging is corrugated paper.)

Warning: Don't place plastic bags within the reach of children, so as to prevent suffocation accidents.

- If it is necessary to move the freezer before it is unpacked, use a forklift or a clamping cart. If a forklift is used, insert the forklift through the bottom of the wooden pallet on the front or back of the case to move the case.
 - If a gripper truck is used, it is clamped from the bottom of the wooden pallet and only the sides of the case are clamped.
 - If the freezer is to be moved after the packaging has been removed, it may be pushed with the help of the equipment's casters.
 - Do not tilt the freezer at an angle greater than 45° during handling.
2. Inventory of accessories: Check the accessories and materials according to the packing list.
 3. Cleaning: Clean the product once before use.
 4. Adjust anti-skid leveling feet: Use a wrench to rotate the fixing foot clockwise so that the fixing foot extends outward and is supported on the ground. Ensure that the equipment will not move when in use.

4.5 First Power-on

When using the equipment for the first time, follow these steps:

1. After the equipment is placed, leveled and cleaned, let it stand for more than 24 hours, and then power it on to ensure the normal operation of the equipment.
2. Under no-load condition, connect the power cord to a special outlet with appropriate specifications.
3. After switching on, turn on the power switch of the equipment (if the auxiliary cooling system is installed, turn off the switch first).
4. Set the required temperature in the freezer (in order to give full play to the best performance of the freezer, if the user needs to use the freezer at a storage temperature below -80°C, set the temperature at -80°C when using it for the first time. After the system runs stably at -80°C for 3~4 hours, the temperature can be further lowered to below -80°C).
5. In case of an alarm after power-on for the first time, click the mute button to eliminate the alarm sound.
6. Check whether the operating temperature of the equipment reaches the required value, observe the normal start and stop of the equipment for more than 24 hours, and put a small amount of items in the freezer after confirming the normal performance.
7. Store items in batches, with the items not exceeding 1/3 of the freezer volume each time. Ensure that the equipment is running properly after shutdown for more than 12 hours before putting in the next batch of items).
8. Try not to open the door during cooling, otherwise the temperature will rise.
9. Turn on the auxiliary cooling system if it is installed.

4.6 Operation after Power Failure

1. The equipment has a memory function for the set value. When the power is restored after power failure, the equipment will continue to operate before power failure.
2. Once the equipment is powered off, it requires 5 minutes before it can be switched on again, so as to avoid damaging the compressor.
3. We guarantee the normal operation of this equipment under certain conditions, but we are not responsible for any loss or damage of stored items after power failure.

Notes:

- The equipment is designed for transportation and relocation using its integrated casters.
- A designated person should be responsible for daily monitoring and recording of the equipment's operating status (recommended inspection every 2–4 hours). In the event of a malfunction or shutdown, the internal temperature will rise. If the unit cannot be repaired promptly, remove the stored items and transfer them to a storage environment that meets their required temperature conditions to prevent damage.
- Before placing items into the equipment, confirm that the operating temperature range of the unit meets the temperature requirements of the stored materials. Failure to do so may result in damage to stored items. When loading the unit, ensure that the air inlet and air outlet are not blocked.
- Due to refrigeration system inertia, a difference may exist between the displayed temperature and the set temperature. This is a normal operating condition.
- This equipment is intended for storage only and must not be used for routine production operations. Do not load large quantities of warm items at one time. Excessive thermal loading may cause prolonged compressor operation and overheating, potentially resulting in compressor damage. Items must be loaded in batches to allow the freezer to cool gradually until the required storage temperature is reached.
- Do not use electrical appliances inside the equipment unless they are approved and certified for such use.
- During operation, short-term power supply fluctuations or brief alarm conditions will not affect normal freezer operation if the voltage variation remains within $\pm 10\%$ of the rated value.
- Do not change the set temperature frequently within short time intervals, as temperature inertia may prevent the system from achieving the desired setpoint. When loading items, ensure adequate air circulation space around the cabinet interior, and do not block the internal temperature sensor, as this will affect temperature stability and control accuracy.
- When loading items, ensure they are properly sealed. Excessive or insufficient moisture content in stored items may affect humidity levels inside the cabinet. Ambient humidity and environmental conditions also influence internal humidity, particularly if the door is opened frequently or not fully closed.


Warning: Children must not play with this equipment. Failure to comply may result in injury or damage.

4.7 Air Switch

Air switch can not only automatically disconnect the power circuit when a machine malfunctions, but also in the event of a short-circuit fault or severe overload to protect the equipment and the power circuit.

5. Location of Parts

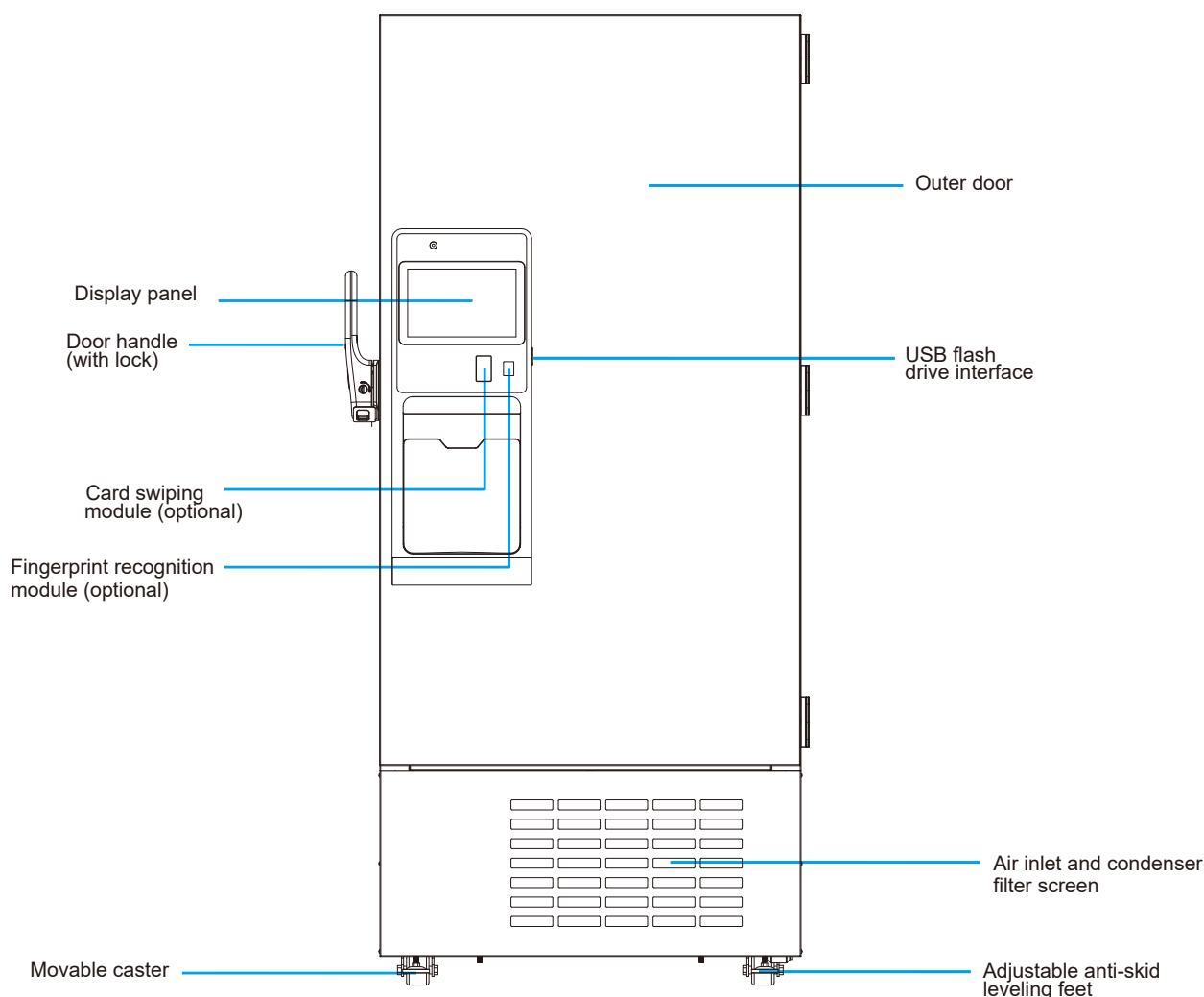* Due to ongoing product improvements and model variations, the actual product may differ from the illustrations shown. Refer to the actual product for details. Illustrations are provided for functional reference only.

Structure and Composition:

The product consists of a freezer cabinet, door (glass door or insulated foam door), refrigeration system, and control system.

Intended Use:

The ultra-low temperature freezer is intended for the cryopreservation and storage of plasma, cells, tissues, and similar materials. It is suitable for use in blood centers, hospitals, epidemic prevention centers, research institutes, and similar facilities.

Intended Users:

Qualified professionals working in hospitals, epidemic prevention centers, research institutes, or similar professional environments.

Medical Conditions:


The ultra-low temperature freezer is intended solely for the storage of plasma, cells, tissues, and similar materials. It is not intended to diagnose, treat, or monitor any medical condition.

Indications for Use:

The ultra-low temperature freezer is intended for the storage of plasma, cells, tissues, and similar materials. There are no specific clinical indications associated with its use.


Contraindications:

None.

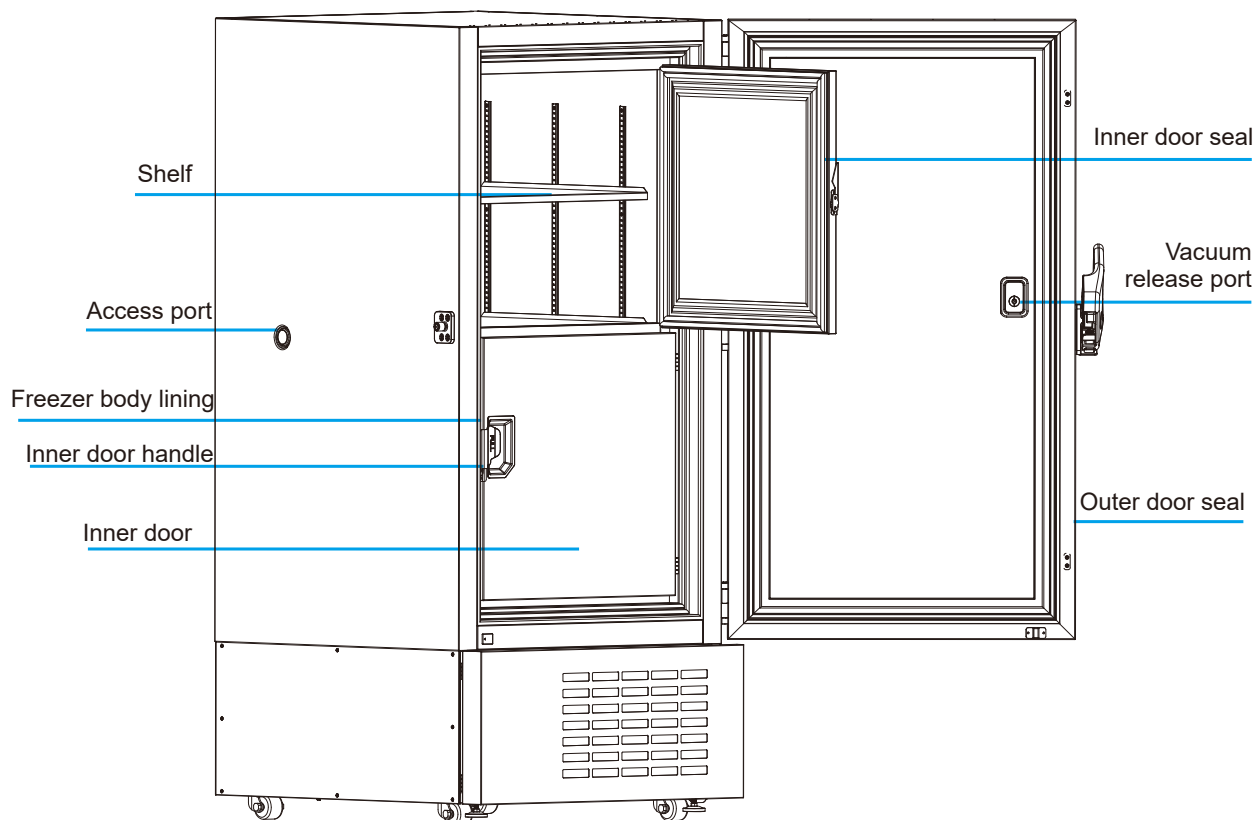! **Note:** For products equipped with the safety function of the electric control lock, log in to release the electric control lock before opening the door.

Open the
filter screen
from here

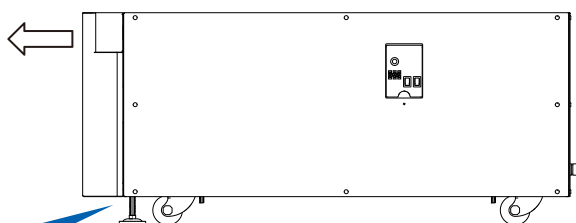
Adjustable anti-skid
leveling feet

- There are adjustable anti-skid leveling feet in front of the two front casters at the bottom of the equipment. Please use accompanying tools  to turn the leveling feet downwards after the equipment is placed securely, so as to stabilize the machine.
- When you need to move the equipment, turn the leveling feet upwards. Do not move it forcibly so as to avoid damage!


6. Operating Instructions

6.1 Introduction to the Homepage:

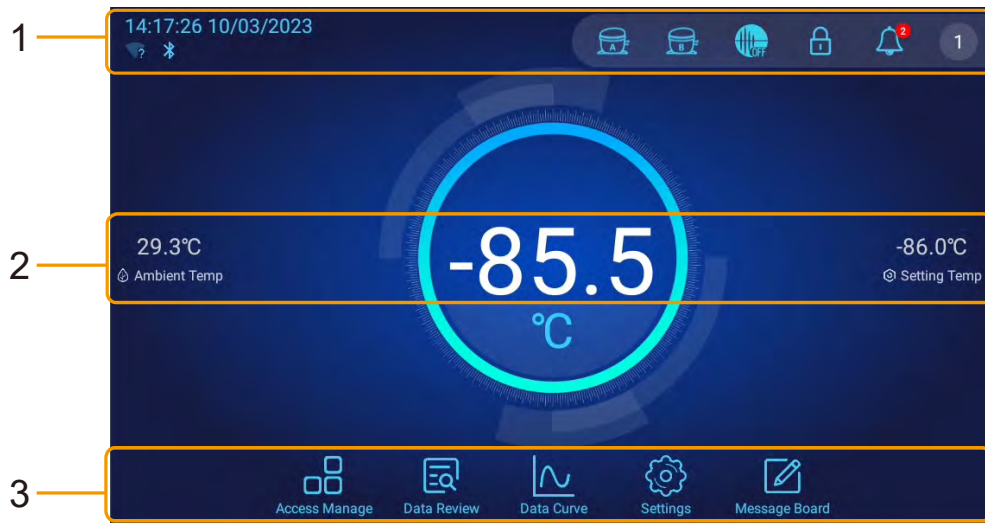
1. Status Bar

a. Time, Date, and Connectivity

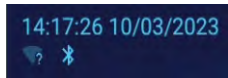The upper left section displays the current time and date as well as the Wi-Fi and Bluetooth connectivity status.

b. Compressor Status Indicators

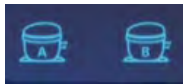- Blue: Normal compressor operation
- Yellow: Compressor stopped
- Red: An error was detected on one or more compressors

c. Mute

- ON – Alarm is muted
- OFF – Alarm is audible

Note: The system may automatically exit mute mode if required to maintain freezer performance or operational stability.

d. Electronic Lock

Tap this icon to unlock the electronic door lock. The electronic lock must be deactivated before opening the freezer door.

- Unlocking the door requires user login. If the user is not logged in, the door will not unlock.
- The door will re-lock automatically and does not require manual locking.
- If the door remains open for an extended time, the handle will lock and prevent the door from closing. To unlock, tap the electronic lock icon again then close the door.

e. Alarm Status

- When an alarm activates, an audible alarm will sound, and a red status indicator will flash. The flashing will continue till the fault is resolved.
- Tap the bell icon once to turn the alarm ON or OFF.
- Press and hold the icon to access the alarm management interface (login is not required).
- If there are unread alarm messages, the number of unread messages is displayed in the upper-right corner of the icon.

2. Temperature Information Display


a. Current Internal Temperature

- When the internal temperature is within the temperature range, the outer ring displays blue-green gradient.
- If the temperature is outside the alarm range, the outer ring displays red-orange gradient.
- Tap the temperature value to enter the operating status interface (login not required).

b. Ambient Temperature

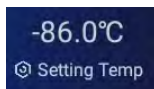
c. Temperature Setting

3. Function Bar

Access Management:

To access and manage stored items. Login is required.

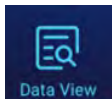Data View:

To display operating history and data usage. Login is not required.

Data Curve:

To view temperature trend curves over time. Login is not required.

Settings:

To configure parameters and functions. Login is required.


Message Board:

To view and post messages with other users.

The number of unread messages will appear on the upper-right corner of the icon. Login is required.

6.2 Alarm Management

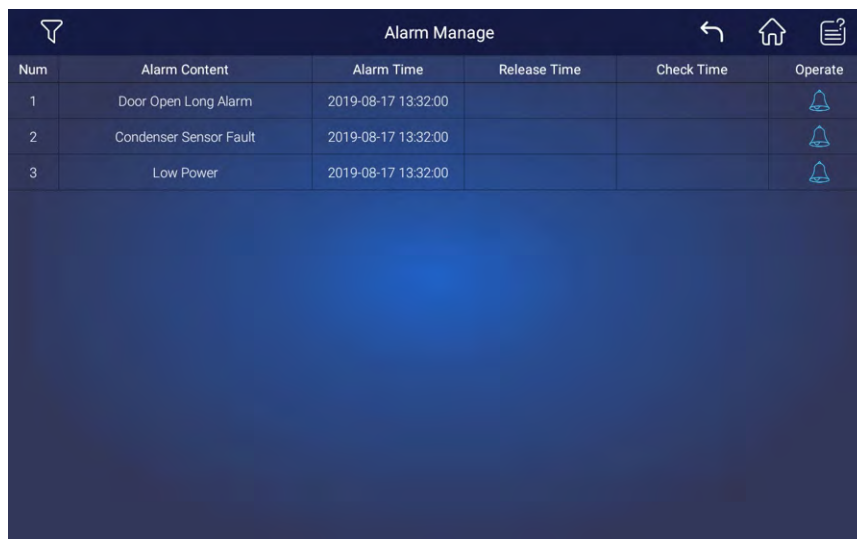The screenshot shows the 'Alarm Manage' interface. At the top, there is a title bar with a filter icon on the left and back, home, and help icons on the right. Below the title bar is a table with the following data:

Num	Alarm Content	Alarm Time	Release Time	Check Time	Operate
1	Door Open Long Alarm	2019-08-17 13:32:00			
2	Condenser Sensor Fault	2019-08-17 13:32:00			
3	Low Power	2019-08-17 13:32:00			

The Alarm Management displays existing and previous alarm information. Users may review, acknowledge, and filter alarm records.

The alarm list includes the following information:

- **Alarm Content:** Description of the alarm event
- **Alarm Time:** Date and time the alarm occurred
- **Release Time:** Date and time the alarm was cleared
- **Check Time:** Date and time the alarm was acknowledged

Tap the filter icon in the upper-left corner to filter alarm records by date.

Icon Descriptions

Acknowledge Alarm

Marks the selected alarm entry as viewed.

Back

Returns to the previous screen.

Home

Return to the main interface.

Help

Opens the help information for the Alarm Management interface.


The alarm functions of this equipment are as follows:

Alarm items	Alarm description
Communication Fault	When the transmission signal between the display screen and the main control board is interrupted, a communication failure alarm will activate.
Door Opening Timeout	When the door has been open for more than 1 minute
Ambient Sensor Fault	When the ambient temperature sensor fails
Power Off Alarm	When the main power supply fails
Low Battery Alarm	When the battery is too low
Temperature Sensor Fault	When the temperature sensor fails
High Temp Alarm	When the temperature in the cabinet is higher than the set temperature of high temperature alarm
Low Temp Alarm	When the temperature in the cabinet is lower than the set temperature of low temperature alarm
Ambient High Temp Alarm	When the ambient temperature is too high
Battery Power Check Fault	Battery power detection failure
Pressure High Alarm	When the pressure of the refrigeration system is too high

6.3 Access Management

a. Take Out Management:

1. Tap **Take Out**
2. Use + (add) or - (remove) to select the quantity.
3. Tap **Take** on the upper-left corner to confirm.

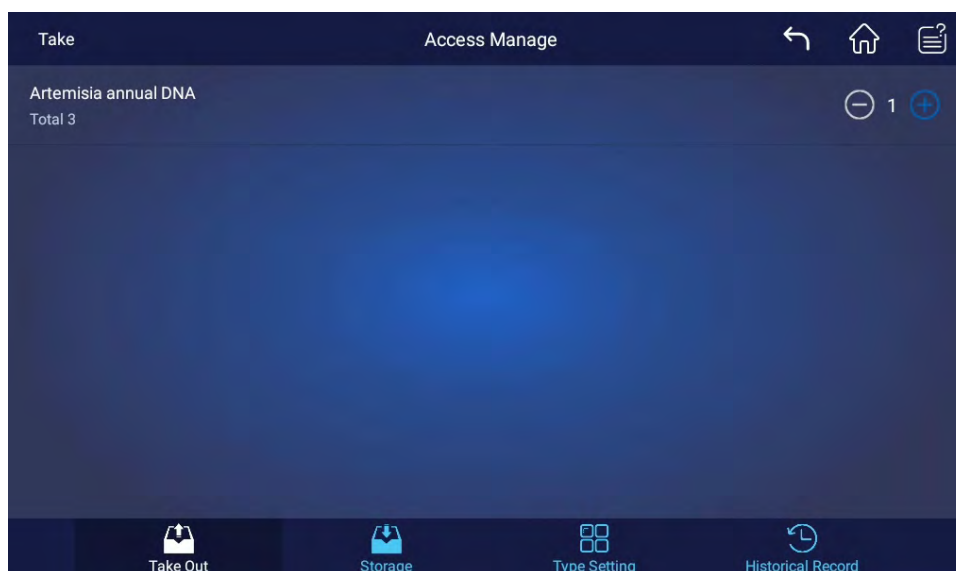b. Storage Management:

1. Tap **Storage**
2. Use the drop down list to select the item type
3. Tap the + (add) or - (remove) to select the quantity
4. Tap **CONFIRM** to adjust

Access Manage

Select Type: Artemisia annual DNA

Select Quantity: 1

The current inventory is 3

CONFIRM

Take Out Storage Type Setting Historical Record

c. Type Setting

1. Tap **Type Setting**
2. Tap **Add** on the upper-left corner to create a new type
3. To edit an existing type, tap **Operate** to the right of the corresponding item type

Add


Access Manage

Num	Type Name	Status	Operate
1	Artemisia annual DNA	Enable	Operate

Take Out Storage Type Setting Historical Record

d. Historical Record

1. Tap **Historical Record**
2. Tap the filter icon on the upper left corner to filter records by date.

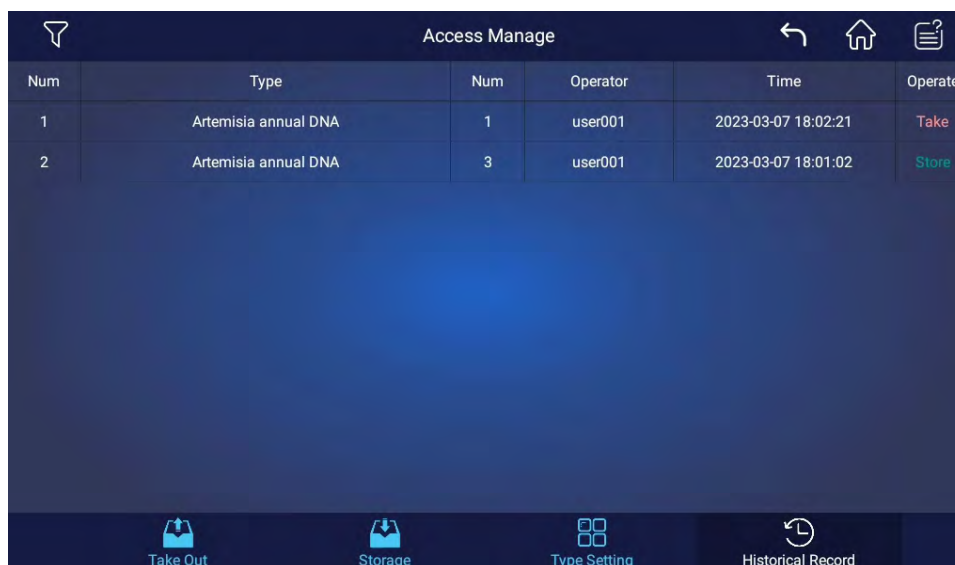Num	Type	Num	Operator	Time	Operate
1	Artemisia annual DNA	1	user001	2023-03-07 18:02:21	Take
2	Artemisia annual DNA	3	user001	2023-03-07 18:01:02	Store

Bottom navigation bar: Take Out, Storage, Type Setting, Historical Record

6.4 Data View

a. Temperature Record

1. Tap Temp Record
2. Tap the filter button in the upper left corner to filter records by date

Num	Time	Internal Temp	Condenser Temp	Condenser2 Temp	Ambient Temp	Door
1	2023-03-10 17:24:00	-85.5℃	29.0℃	28.3℃	28.4℃	Off
2	2023-03-10 17:19:00	-85.5℃	29.0℃	28.3℃	28.5℃	Off
3	2023-03-10 17:14:00	-85.5℃	29.0℃	28.3℃	28.5℃	Off
4	2023-03-10 17:09:00	-85.5℃	29.0℃	28.3℃	28.4℃	Off
5	2023-03-10 17:04:00	-85.5℃	29.0℃	28.3℃	28.4℃	Off
6	2023-03-10 16:59:00	-85.5℃	29.0℃	28.3℃	28.4℃	Off
7	2023-03-10 16:54:00	-85.5℃	29.0℃	28.5℃	28.5℃	Off
8	2023-03-10 16:49:00	-85.5℃	29.0℃	28.5℃	28.6℃	Off
9	2023-03-10 16:44:00	-85.5℃	29.0℃	28.5℃	28.5℃	Off
10	2023-03-10 16:39:00	-85.5℃	28.9℃	28.6℃	28.7℃	Off
11	2023-03-10 16:34:00	-85.5℃	28.9℃	28.7℃	28.7℃	Off
12	2023-03-10 16:29:00	-85.5℃	28.9℃	28.7℃	28.7℃	Off
13	2023-03-10 16:24:00	-85.5℃	28.9℃	28.7℃	28.6℃	Off
14	2023-03-10 16:19:00	-85.5℃	28.9℃	28.8℃	28.5℃	Off

Bottom navigation bar: Temp Record, Operation Record, Alarm Record, Detailed Record

b. Operation Records

- 1. Tap Operation Records
- 2. Tap the filter button to filter records by date

Note: The operation record types includes: login, logout, user addition, user modification, export, temperature setting, alarm setting, unlock, printing, door opening, and door closing.

Data View

Num	Time	Operator	Operation Item(Par)	Operation Item(Sub)
1	2023-03-10 17:24:57	user001	Authority	Unlock
2	2023-03-10 17:24:47	user001	Authority	Login
3	2023-03-10 14:28:02	user001	Authority	Logout
4	2023-03-10 14:06:16	user001	Authority	Login

Temp Record

Operation Record

Alarm Record

Detailed Record

c. Alarm Record

Tap alarm record to view the alarm history

Data View

Num	Alarm Content	Alarm Time	Release Time
1	High Temp Alarm	2023-03-10 14:09:00	
2	Battery Powered Alarm	2023-03-10 14:09:00	

Temp Record

Operation Record

Alarm Record

Detailed Record

To Export Data on to a USB

1. Tap the Export button to open the export range selection dialog
2. Select the start date and end date
3. Tap **CONFIRM** to export the data as an Excel file (*.xls) to the USB flash drive


To Print Data

Tap the Printer button to open the print settings dialog, where you can choose Print By Time or Print By Day.

- For Print by Time: Set the start time, end time, and print interval.
- For Print by Day: Set the start time, end time, and the specific time of day for printing.

Note: A printer must be connected to use this function. The equipment supports connection to a Bluetooth printer via Bluetooth.

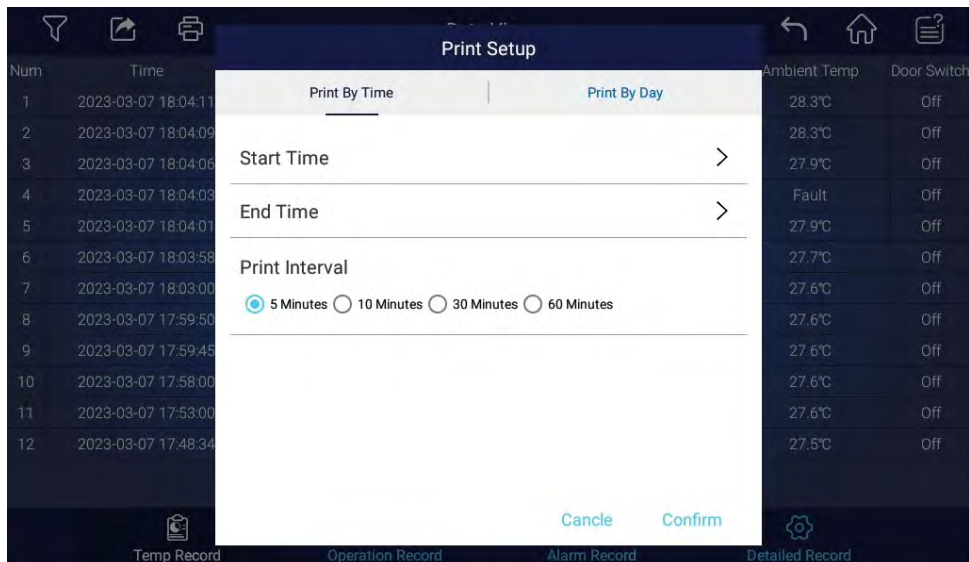6.5 Data Curve


- Tap on the **Day/Week/Month** tab to view the curve graph of the corresponding tab
- Tap the **Previous Day/Next Day** tab to switch to the curve graph corresponding to the date
- Tap the filter button in the upper-left corner to jump to the curve graph of the corresponding date
- Tap the export button in the upper-left corner to export the curve figure to the USB flash drive
- Pinch-to-zoom to zoom in/out the curve graph, and slide with a single finger to view detailed temperature value

6.6 Settings


a. Temperature Setting

Adjust the set temperature by sliding the temperature control or tapping the temperature value directly. Tap °C (Celsius) or °F (Fahrenheit) to switch the temperature unit.

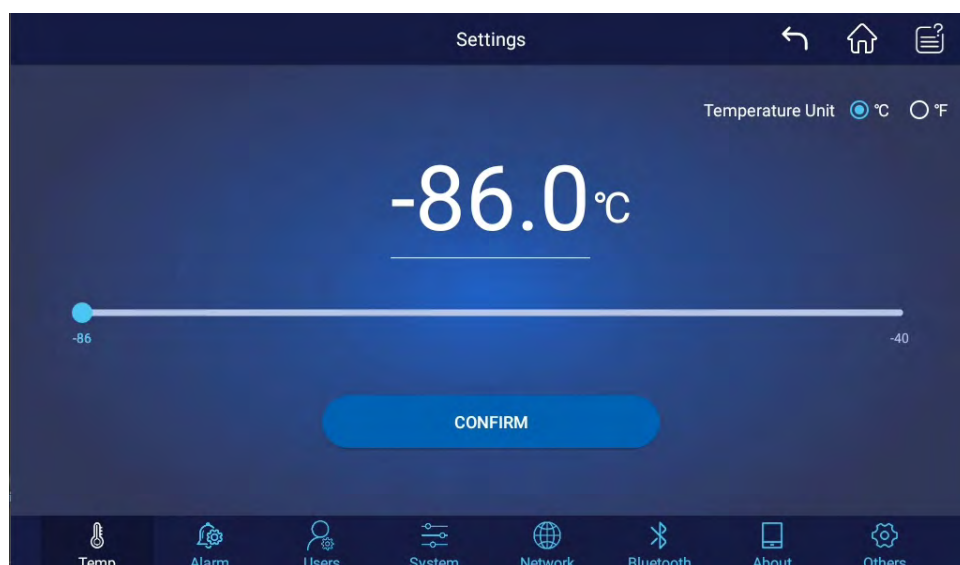
b. Alarm Setting

To adjust the High-Temperature Alarm , Low-Temperature Alarm , and the High Ambient Temperature Alarm values, use the slider or tap the values to adjust.

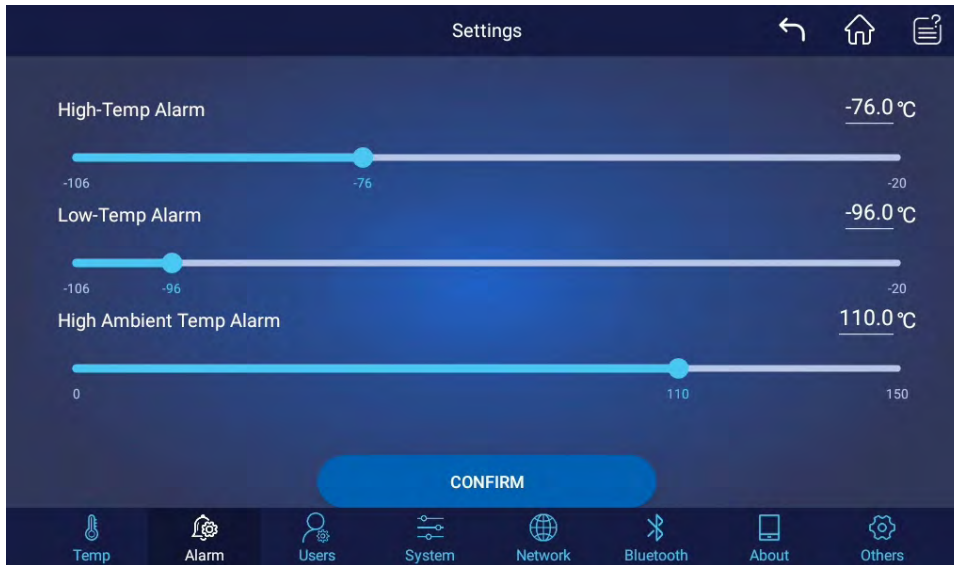
c. User Management Settings

- The administrator can add or modify users.
- The user created by the administrator are ordinary user and do not have the authority to add/modify contents.
- Tap **Reset Password** to reset the user's login password. The default password is **123456**.

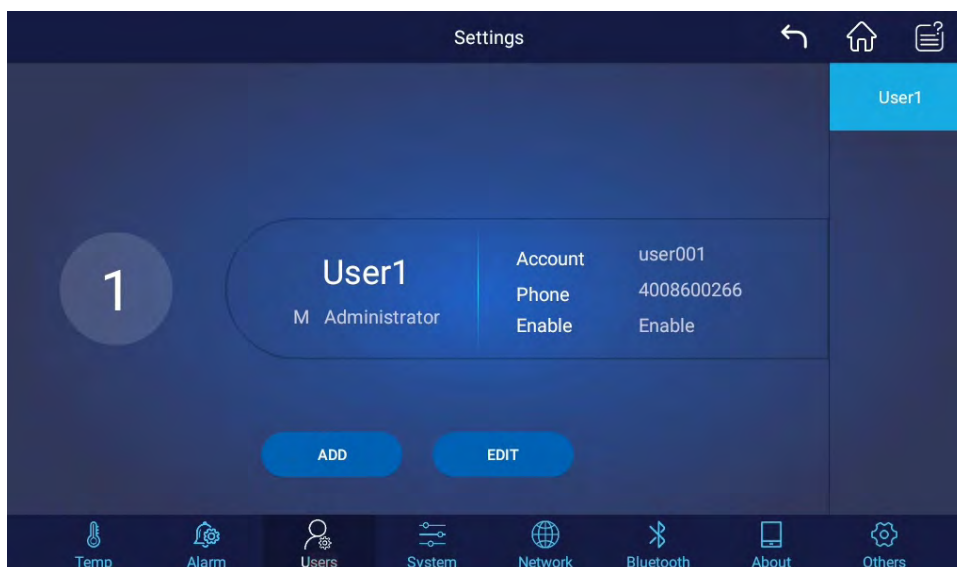
d. System Settings

The System Settings allow you to adjust the display brightness, screen timeout (automatic screen off when there is no operation), the time for the display to automatically return to the main interface when there is no operation, volume level, alarm sound (changes take effect after returning to the main interface), date/time, automatic user logout, and language.

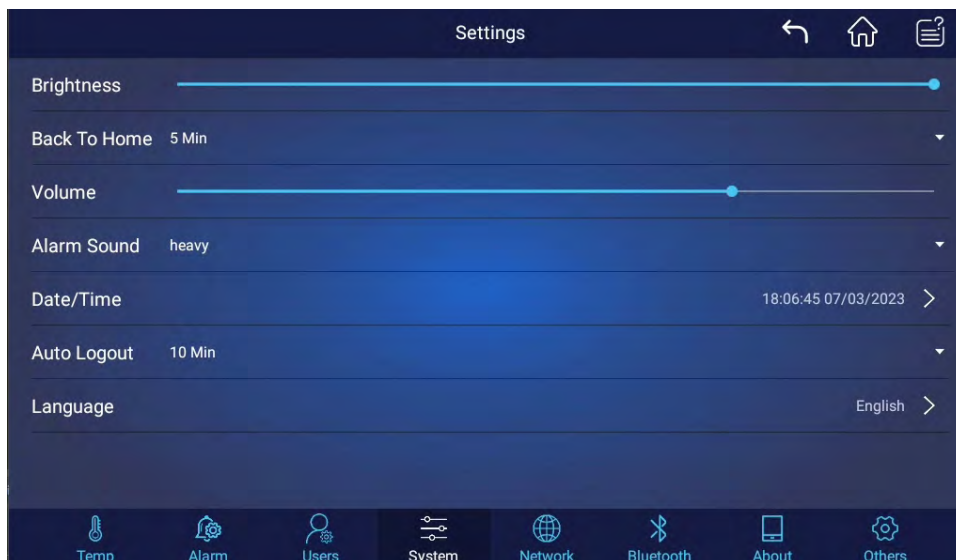
e. Network Settings

Set up the Wi-Fi connection to link to a wireless network.

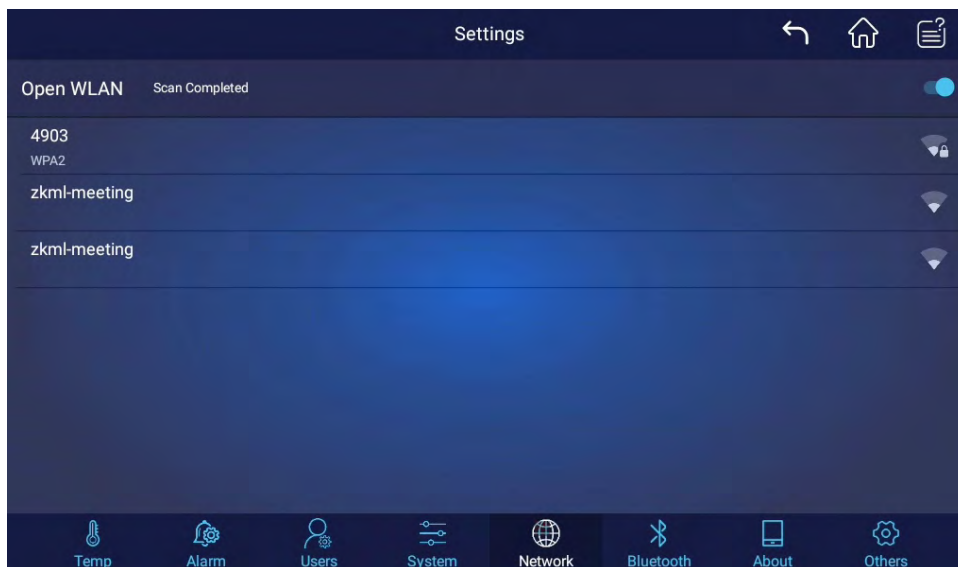
f. Bluetooth Settings

You can set up the Bluetooth connection.

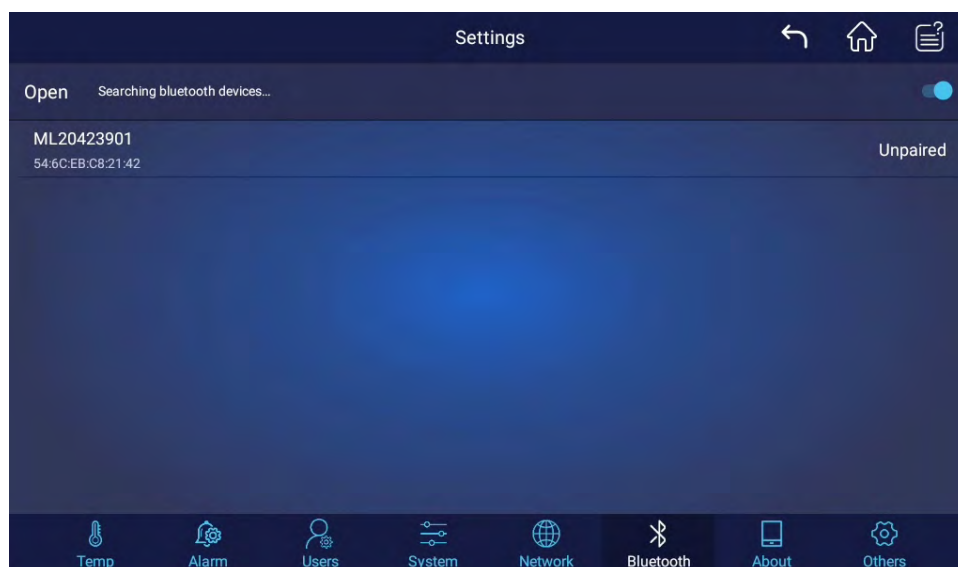
g. About Screen

- This page displays the current machine model, software version, and hardware version.
- If the machine is connected to the Internet, tap **Check Update** to check for available software updates.
- When a newer version is available, you can choose to download and install the update.

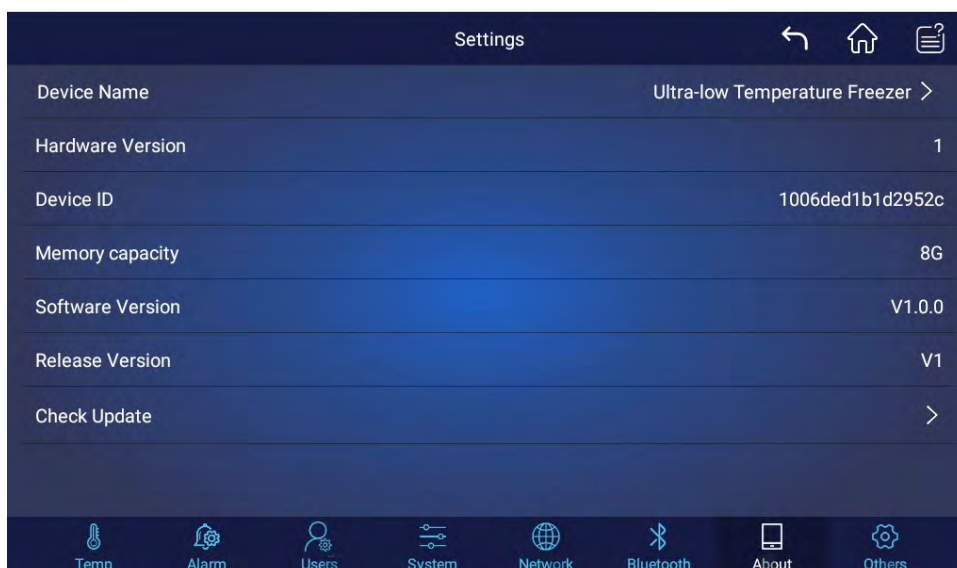h. Other Settings

This page is accessible only to administrators. From this page, the administrator can configure the display of virtual buttons (show or hide), set the initial password for new user logins, enable or disable encryption for exported Excel files and define the encryption password, set the allowable temperature range, select the printer, and configure the equipment name.

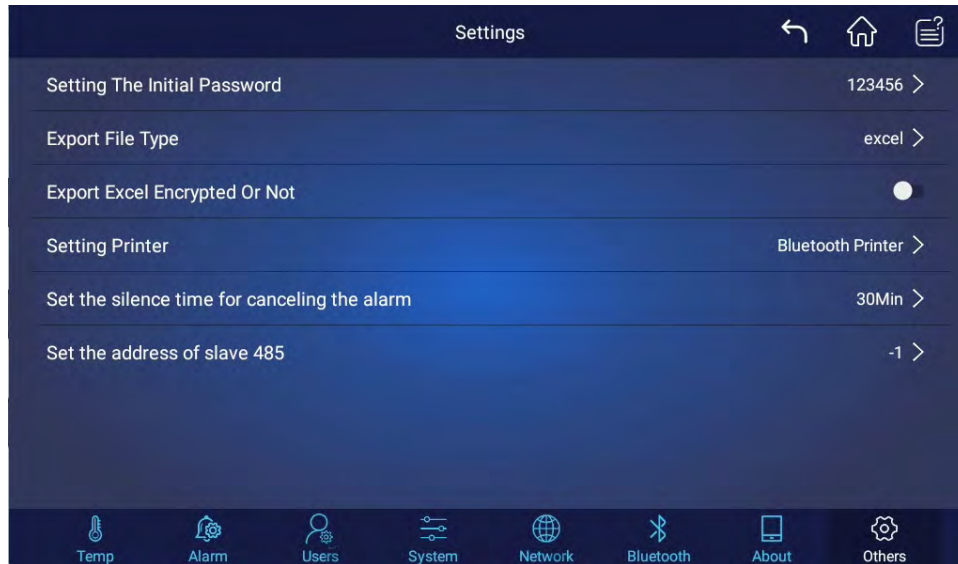6.7 Message Board


- Tap **Add (+)** on the lower-right corner to open the message posting dialog.
- You can choose to make the message visible to All Users or Specific Users.
- Selecting All Users allows all users to view the message, while Specific Users restricts visibility to the selected users.

Tap the **Message** icon to mark the message entry as viewed.


Tap **Mine / All** on the upper-left corner to switch between your messages and all messages.

To delete the message, tap and hold the message entry and slide to the left. The delete button will appear on the right side of the message.

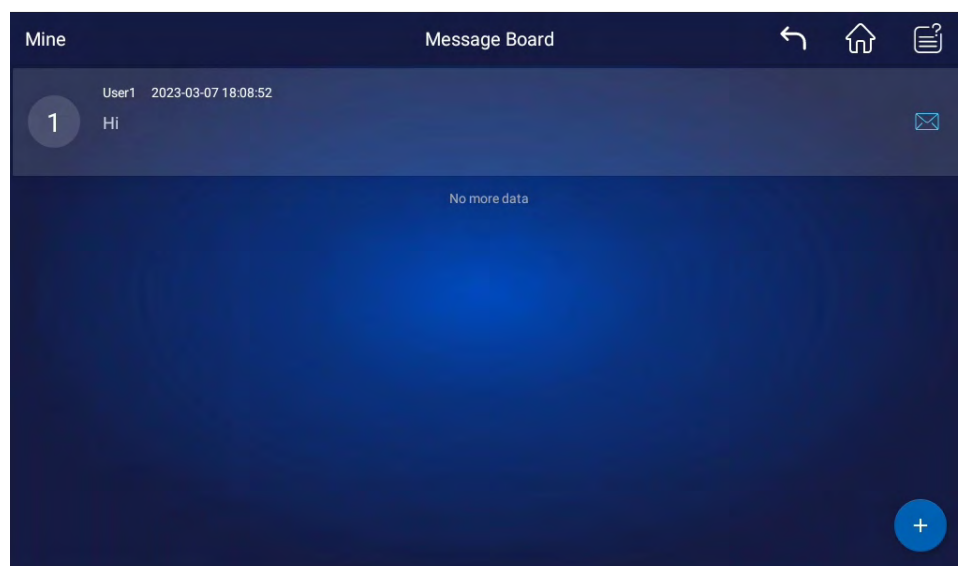6.8 Login Options

a. Account and Password

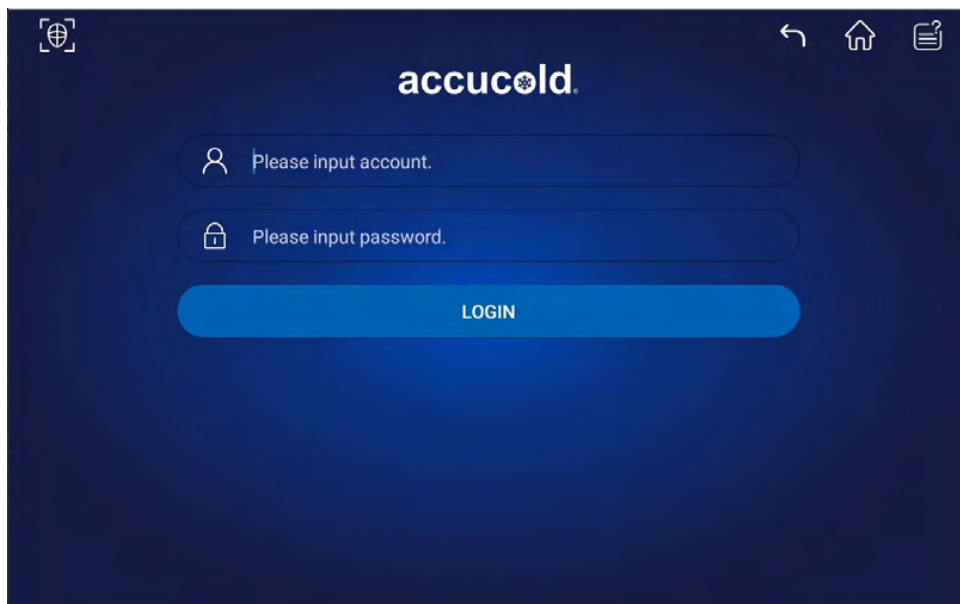The image shows the login screen of the 'accucold' application. At the top, there is a navigation bar with a globe icon on the left and back, home, and help icons on the right. The title 'accucold.' is centered. Below the title, there are two input fields: the first is labeled 'Please input account.' with a person icon, and the second is labeled 'Please input password.' with a lock icon. A blue 'LOGIN' button is positioned below these fields.


Resetting the password:

Enter the correct account number and password to log in. (The default password for new users is **123456**).

You will be asked to change the password when logging in for the first time.

The administrator can change the default password for new users at anytime.


1. Tap Settings
2. Other settings
3. Set the initial password

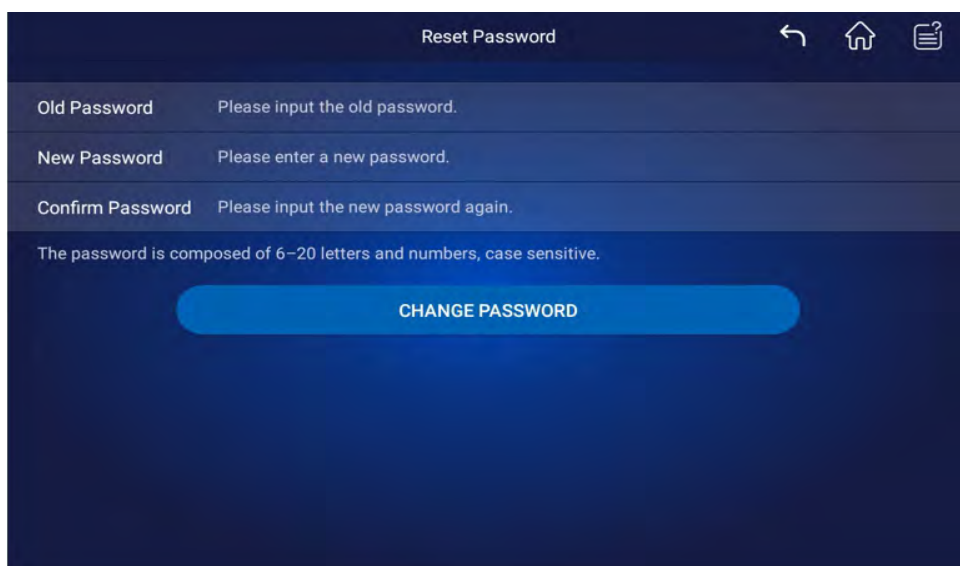The image shows the 'Reset Password' screen. At the top, there is a navigation bar with the title 'Reset Password' and back, home, and help icons. Below the title, there are three input fields: 'Old Password' with the prompt 'Please input the old password.', 'New Password' with the prompt 'Please enter a new password.', and 'Confirm Password' with the prompt 'Please input the new password again.'. Below these fields, there is a note: 'The password is composed of 6-20 letters and numbers, case sensitive.'. A blue 'CHANGE PASSWORD' button is at the bottom.


b. Fingerprint (optional)

1. Login as the user
2. Select the Fingerprint Registration option on the homepage menu (top right)
3. Register the fingerprint

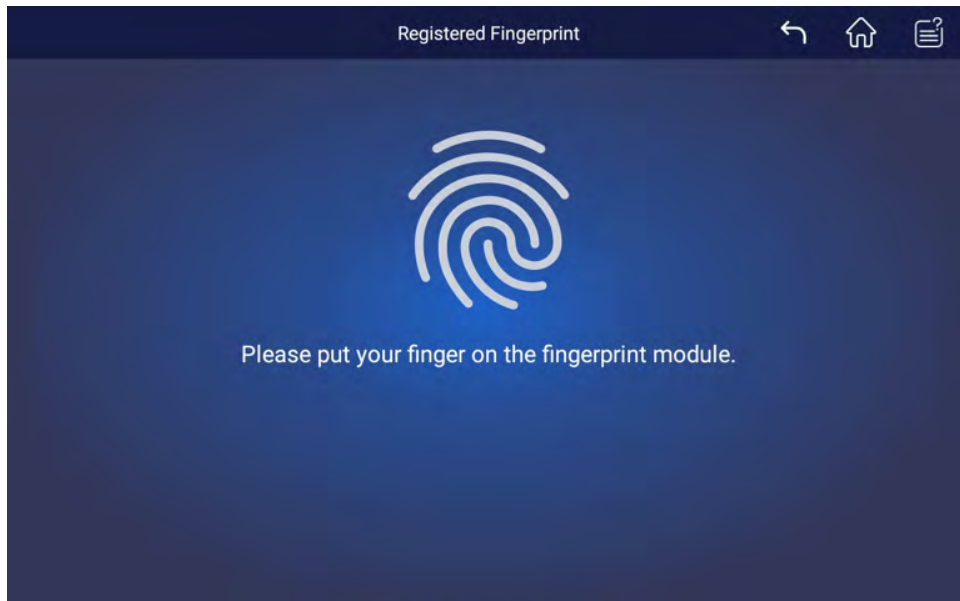c. Identification Card (optional)

1. Login as the user
2. Select the Fingerprint Registration option on the homepage menu (top right)
3. Register the fingerprint

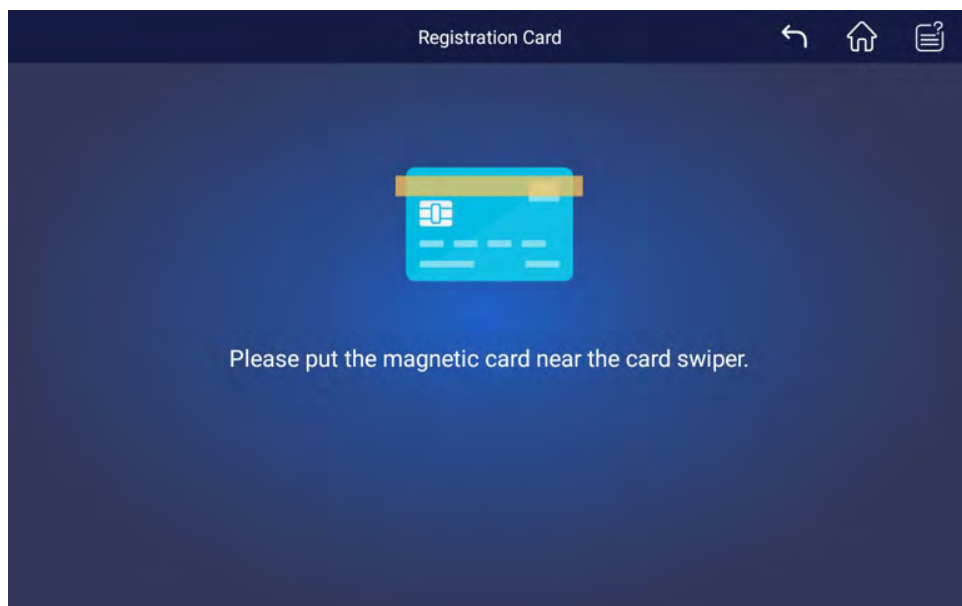d. Facial Recognition (optional)

To Register Face ID:

1. Login as the user
2. Select the Facial Registration option on the homepage menu (top right)
3. Once recognized, tap the  icon to register the face
4. To use the facial recognition, tap the  icon on the top left of the login screen

Note: To use facial recognition offline, the face ID must be registered when the unit is online. If the face cannot be registered after the machine is powered off, it can be activated again online.

6.9 System Features

a. Freezer Rack and Freezer Box


Model	Internal Dimensions of the Freezer Body (W x D x H)	Shelf/Inner Door	Freeze Racks	Dimensions of the Freeze Rack (W x D x H)	Freezer Box Max. Capacity
AULT530DC	585 x 696 x 1266 23" x 27 ³ / ₈ " x 49 ⁷ / ₈ "	3/2	4×4 pcs	140 x 685 x 293 5 ¹ / ₂ " x 27" x 11 ¹ / ₂ "	400 pcs 5 x 5 x 16 cm 2" X 2" X 6 ¹ / ₄ "
AULT780DC	865 x 696 x 1286 34" x 27 ³ / ₈ " x 50 ⁵ / ₈ "	3/2	4×6 pcs	140 x 685 x 293 5 ¹ / ₂ " x 27" x 11 ¹ / ₂ "	600 pcs 5 x 5 x 24 cm 2" X 2" X 9 ¹ / ₂ "

b. Code Scanning Window


In any operation interface, you can scan a QR code or a bar code in the built-in code scanning window of the product to fill the information encoded in the QR code or the bar code in the current cursor position.

6.10 Remote Alarm Terminal and RS-485 Interface

The remote alarm terminal is installed on the back or the lower part of the right side of the freezer body, and the alarm signal is output by the terminal. The load-bearing capacity of the terminal is **30V, 2A DC**. Users can choose normally open (N.O.), normally closed (N.C.) and common terminals (COM) according to their own needs.

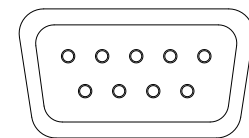RS-485

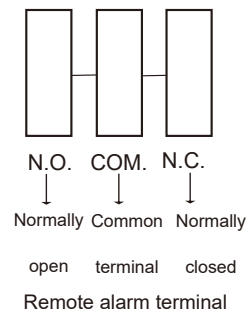7. Defrost, Maintenance, and Discontinuation

7.1 Defrost and Cleaning

Defrosting:

Cooling efficiency reduces and power consumption increases when there's excessive frost buildup on the inner surfaces of the unit. Use the defrosting tool supplied with your unit to remove ice as needed. Do not use sharp or metal tools, as they may damage the equipment.

Before defrosting:

- Remove all stored items
- Disconnect the power supply

After defrosting is complete:

- Reconnect the power
- Allow the freezer to reach the set temperature
- Returning items into the cabinet

Always wear appropriate protective gloves when loading or unloading items to prevent frostbite.

Cleaning and Maintenance:


- Clean and maintain the equipment regularly.
- For safety, unplug the power cord before cleaning.
- Wipe both the interior and exterior surfaces of the freezer using a soft cloth dampened with warm water.

Cleaning the Filter Screen:

Remove dust from the filter screen regularly to ensure proper heat dissipation.

Cleaning Procedure:

1. Pull the upper-right corner of the front ventilation panel outward to open the panel.
2. Lift and remove the filter screen.
3. Remove surface dust using a soft brush or rinse with tap water.
4. Allow the filter screen to dry completely before reinstalling it.

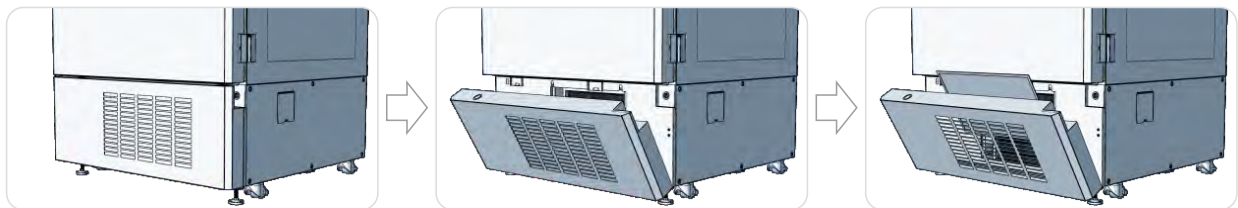Warning:

- To prevent electric shock or personal injury, always disconnect the power supply before performing any repair or maintenance.
- During maintenance, avoid inhaling chemicals, drugs, or airborne particles in the surrounding area, as inhalation may be harmful to your health.

Notes:

- Do not spray or pour water directly onto the cabinet, as this may reduce the insulation performance of electrical components and cause corrosion of metal parts.
- Do not use hot water, corrosive cleaners, or organic solvents to clean the freezer.
- Do not place heavy objects on top of the equipment, as excessive weight may cause deformation or damage.
- The air circuit breaker supplied with the equipment supports both manual operation and automatic protection. In the event of an electrical short circuit or overload, the circuit breaker will automatically cut off the power supply to protect the equipment and the electrical system.

7.2 Equipment Discontinuation

Deactivation:

The unit is sitting idle for an extended period:

- Disconnect the power supply
- Clean the interior and exterior surfaces with a soft cloth dampened with warm water
- Allow the unit to dry thoroughly and seal it
- The equipment must be locked to prevent children from opening the freezer door

Disposal:

When the equipment reaches the end of its service life, it must be disposed of in accordance with local regulations and transferred to a qualified recycling or disposal agency. Unauthorized personnel must not disassemble or dismantle the equipment. Discarded equipment must be placed in a designated area inaccessible to children to prevent injury or other hazards.

7.3 Battery Maintenance:

1. To extend battery service life and prevent degradation due to prolonged inactivity, operate the equipment for at least 24 hours each month to allow the battery to recharge fully.
2. If the equipment will not be used for an extended period, connect it to the power supply periodically (generally once per month) and turn on the power lock switch so the unit can operate and complete a full charging cycle. Each charging period should last no less than 24 hours.
3. If the main power supply is interrupted, turn off the equipment's power lock switch promptly. Failure to do so may result in battery discharge and, in severe cases, permanent battery damage.
4. During long-term power outages or transportation, the main power switch must be turned off. Otherwise, prolonged battery discharge may occur, potentially leading to permanent battery damage and abnormal display behavior when power is restored.
5. The battery is a consumable component with an expected service life of approximately 2–3 years. Improper use — such as repeated deep discharge or end-of-life degradation — may trigger a low battery alarm. Refrigeration performance is not affected; however, alarm functions, including power failure alarms, may not operate correctly.

For additional assistance, contact our Customer Service facility at **800-932-4267** (Ext. **513**) between 9:00 AM and 5:00 PM ET or visit our website: **summitappliance.com/support** at any time for your battery replacement.

Replacing and Recovering Battery

1. Turn off the power switch and unplug the power cord from the electrical outlet.
2. Detach the right-hand panel of the freezer's cabinet frame with a screwdriver.


Note: The battery can be replaced without opening the electrical box. The electrical box contains high-voltage components and must only be opened by qualified engineers or authorized service personnel. Before any access, ensure the power supply is turned off, the power plug is unplugged, and the power switch is in the OFF position.

3. Locate the battery on the left side, outside the electrical box and disconnect the battery cable.

Important: Before disconnecting, ensure to note the polarity and cable connection order to prevent damage to the control system when installing the new battery. Reversed polarity may cause permanent system damage. Typically, the red wire connects to the positive terminal, and the black wire connects to the negative terminal.

4. Use a screwdriver to remove the screws securing the battery mounting plate, then remove the battery.
5. Replace the battery with a **BT-12M4.0AC** (12V, 4.0Ah) battery only.
6. Dispose of the removed battery through a qualified local battery recycling agency in accordance with applicable regulations.

Important: To ensure proper replacement and prevent damage to the control system due to incorrect installation, it is strongly recommended to contact authorized service personnel for battery replacement or technical guidance.

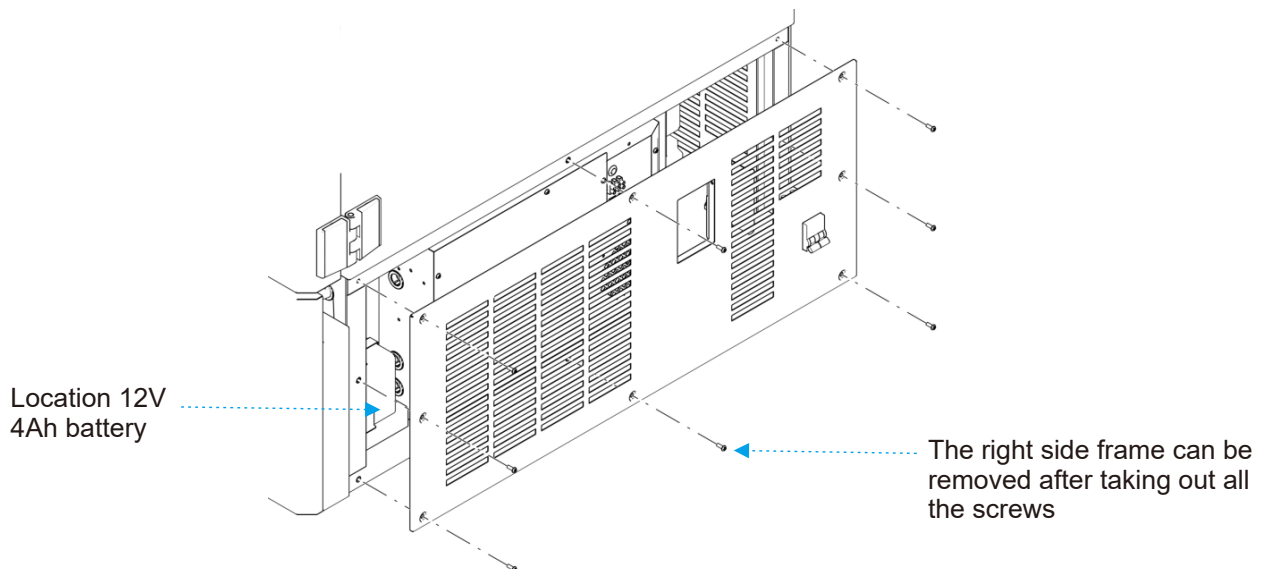Basic Requirements for Operations

- Ensure the operating environment has adequate ventilation.
- Clean the work area to confirm there are no flammable or explosive materials or open flames present.
- Ensure the work site is equipped with a ventilation fan and a dry powder fire extinguisher.
- Prepare a sufficient supply of nitrogen or other inert gas for use.
- Inspect oxygen and acetylene cylinders, pressure regulators, and pipelines to ensure there is no aging, damage, or leakage.
- All operations must be carried out strictly in accordance with approved procedures and performed only by personnel with the required work permits.
- Wear appropriate personal protective equipment (PPE) correctly, and secure cuffs tightly to prevent burns or scalds from flames or hot surfaces. During operation, verify that all pipe connections are properly tightened and free from looseness or leakage.

Before Operations

- Before starting work, wear appropriate personal protective equipment (PPE), such as a protective mask. Secure cuffs tightly and do not wear scarves or loose clothing.
- Open windows, doors, and activate fans to ensure adequate ventilation at the work site.
- Confirm that the work area is free of flammable or explosive materials and open flames.
- Check that welding tools, pipelines, joints, and connections are properly tightened and free from leakage.
- Verify that nitrogen or other inert gas cylinders are filled with the correct gas. Pay special attention when using the re-ignition method to confirm that the gas is not oxygen.
- Inspect the AB refrigerant gas to be used and confirm that there is no leakage.

Analysis of Risk Sources and Hazard Factors

- Poor ventilation at the repair site may cause the concentration of combustible materials to exceed safe limits, leading to fire or explosion hazards.
 - a. Uncertainties in the environment of the repair work site.
 - b. Accidents caused by oxygen, acetylene pipeline bursting.
- Hazards resulting from mixing inert gases with combustion-supporting gases during leak testing or power-on testing.
- Personal injury caused by failure to follow operating procedures, including burns or scalds during welding operations.

After Operations

- Working in a poorly ventilated environment is strictly prohibited.
- Using oxygen or other combustion-supporting gases in place of nitrogen or inert gases is prohibited.
- Using unverified or unknown gases is prohibited.
- Replacing the specified mixed refrigerant with other refrigerants is prohibited.
- Using aged, damaged, or leaking hoses and joints is prohibited.
- Working in environments containing flammable or explosive substances is prohibited.
- Smoking during operations is strictly prohibited.

**Warning:**

This product uses a flammable refrigerant. Unauthorized disassembly or repair—especially of the compressor, refrigeration piping, or related components—is strictly prohibited. In the event of a malfunction, repairs must be performed by our Service Center.

**Notes:**

- The equipment may only be repaired, maintained, or modified by engineers certified by our company to ensure proper operation and compliance with applicable safety standards.
- Before contacting a maintenance engineer, please clean and disinfect the equipment. During the warranty period, the Company will not be responsible for faults or damage caused by improper use.
- Refrigeration system leak testing must be performed by a qualified professional at least once every 12 months.
- If a refrigerant leak is detected, immediately disconnect the power supply and notify the Company's after-sales service personnel or other qualified technicians for repair.

Troubleshooting

You can solve many common problems easily, saving you the cost of a possible service call. Try the suggestions below to see if you can solve the problem before calling the servicer.

Problems	Causes and solutions
Ultra-low temperature freezer does not operate	• Ensure the power outlet is energized. • Confirm the power plug is fully inserted and secure. • Check that the power supply has not been disconnected. • Verify that the supply voltage is within the specified range (not too high or too low).
Compressor is not running	• Ensure the switch located at the lower right rear of the freezer is turned ON • Confirm that the temperature setting is correct.
The temperature does not reach the set value	• Check that frost buildup on the inner wall is not excessive. • Ensure the door is fully closed and avoid frequent opening in a short period. • Do not load a large quantity of items at one time. • Make sure the ambient temperature is within the specified operating range. • Ensure the air filter is clean and not blocked.
Excessive noise	• Ensure the freezer is placed on a level surface. • Make sure the freezer is not in contact with walls or other objects.
High ambient temperature alarm	Cause: Ambient temperature exceeds 32 °C (90 °F) or insufficient clearance around the freezer restricts ventilation. Solution: Maintain a minimum clearance of 30 cm (12 in.) between the freezer and surrounding walls. Ensure vents are not blocked. Clean the air filter regularly. Improve room ventilation or use air conditioning during hot conditions.

The following conditions are normal and do not indicate a malfunction:

- When the compressor starts or stops, you may hear a slight impact or clicking sound. This is normal.
- A gentle flowing or water-like sound may be heard as the refrigerant moves through the piping during operation.
- During rainy or humid seasons, frost or condensation may form on the exterior surface of the freezer. This can be wiped off with a dry cloth.
- Allow at least 5 minutes between door openings. Opening the door too frequently allows warm air to enter, and the rapid temperature drop may create negative pressure, making the door difficult to open.
- The inner edge of the liner opening is a stress-relief area designed for ultra-low temperatures. Minor cracks may appear due to repeated temperature changes; this is normal and does not affect operation. If the LCD screen flashes during printing or after extended use, this is also normal and indicates data reading or screen refresh.
- If the display shows an abnormal reading or alarm during startup, the control panel may not have fully read the data. The display should return to normal within a few seconds. If the issue persists, turn off the power and restart the freezer.

If you have checked the information above and you still need help with your appliance, call our Customer Service facility at **800-932-4267 (Ext. 513)** between 9:00 AM and 5:00 PM ET or visit our website **summitappliance.com/support** at any time. We will do our best to answer your questions.

To order replacement parts, visit our website: **summitapplianceparts.com**

WARNING: This product can expose you to chemicals including Nickel (Metallic) which is known to the State of California to cause cancer.

For more information go to www.P65Warnings.ca.gov

Note: Nickel is a component in all stainless steel and some other metallic compositions.

NOTES

LIMITED WARRANTY

FIVE-YEAR LIMITED WARRANTY

Warranty coverage is available within the 48 contiguous United States, along with select areas of Alaska, Hawaii, and Canada. For locations not covered, consumers can either contact a local provider or coordinate with Felix Storch, Inc. to find nearby service options, which may qualify for partial reimbursement.

The warranty for this appliance is valid for five years from the date of purchase, provided the appliance is operated and maintained according to the instructions included with the product. The warrantor will cover the cost of factory-specified parts and repair labor needed to correct defects in materials or workmanship. Service must be performed by a designated service company.

All parts, except the compressor, are warranted for five years against manufacturing defects. Plastic parts, shelves, and cabinets are produced to commercially acceptable standards but are not covered for damages incurred during handling or for breakage.

12-YEAR COMPRESSOR WARRANTY

1. The compressor is covered for 12 years.
2. Replacement does not include labor.

ITEMS WARRANTOR WILL NOT PAY FOR:

1. Service calls to correct the installation of your appliance, to instruct you how to use your appliance, to replace or repair fuses or to correct wiring or plumbing.
2. Service calls to repair or replace appliance light bulbs or broken shelves. Consumable parts (such as filters) are excluded from warranty coverage.
3. Damage resulting from accident, alteration, misuse, abuse, fire, flood, acts of God, improper installation, installation not in accordance with electrical or plumbing codes, or use of products not approved by warrantor.
4. Replacement parts or repair labor costs for units operated outside the 48 contiguous United States, along with select areas of Alaska, Hawaii, and Canada.
5. Repairs to parts or systems resulting from unauthorized modifications made to the appliance.
6. Expenses for travel and transportation for product service in remote locations.
7. The removal and reinstallation of your appliance if it is installed in an inaccessible location or is not installed in accordance with published installation instructions.

DISCLAIMER OF IMPLIED WARRANTIES; LIMITATION OF REMEDIES

CUSTOMER'S SOLE AND EXCLUSIVE REMEDY UNDER THIS LIMITED WARRANTY SHALL BE PRODUCT REPAIR AS PROVIDED HEREIN. IMPLIED WARRANTIES, INCLUDING WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED TO FIVE YEARS. WARRANTOR SHALL NOT BE LIABLE FOR INCIDENTAL OR CONSEQUENTIAL DAMAGES. SOME STATES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES, OR LIMITATIONS ON THE DURATION OF IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS, SO THESE EXCLUSIONS OR LIMITATIONS MAY NOT APPLY TO YOU. THIS WARRANTY GIVES YOU SPECIFIC LEGAL RIGHTS AND YOU MAY ALSO HAVE OTHER RIGHTS, WHICH VARY FROM STATE TO STATE.

CALIFORNIA CARB/SNAP DISCLOSURE

This product uses eco-friendly hydrocarbon refrigerant and fully complies with California CARB regulations.

However, we are required by California Law to provide the following disclosure statement in every product sold in California.

"This equipment is prohibited from use in California with any refrigerants on the 'List of Prohibited Substances' for that specific end-use, in accordance with California Code of Regulations, title 17, section 95374. This disclosure statement has been reviewed and approved by Felix Storch, Inc. and Felix Storch, Inc. attests, under penalty of perjury, that these statements are true and accurate."

This product does **not** use any refrigerants on the 'List of Prohibited Substances'"

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
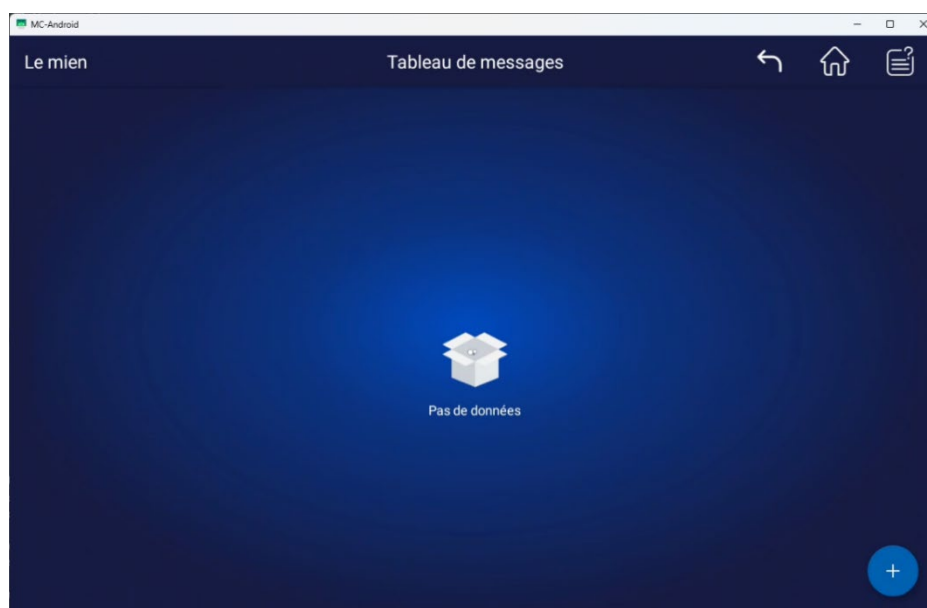
6.8 Options de connexion

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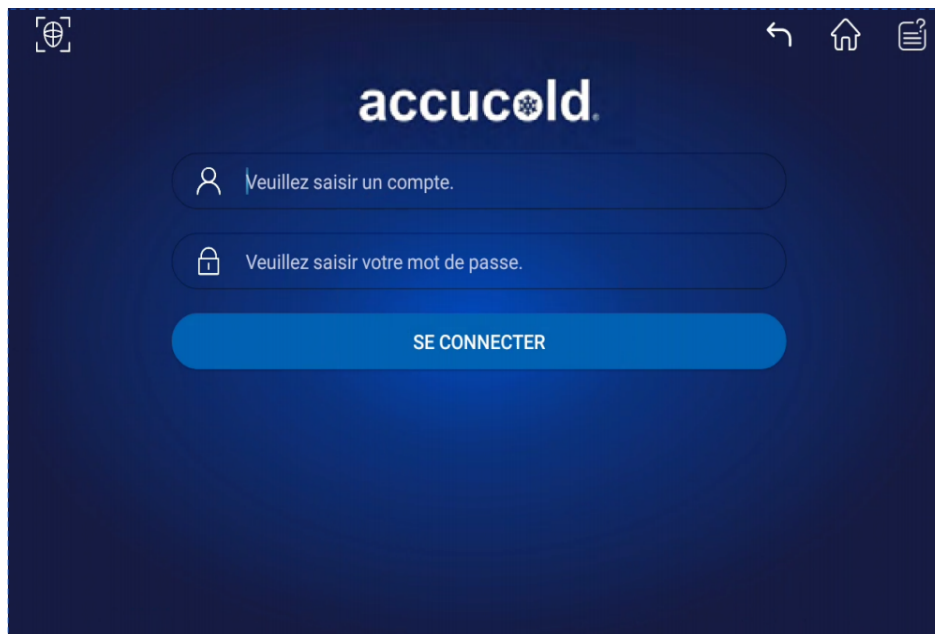The screenshot shows the login interface for 'accucold'. At the top, there is a navigation bar with a globe icon on the left and back, home, and help icons on the right. The 'accucold' logo is centered. Below the logo are two input fields: the first is for the account name with a person icon and the placeholder text 'Veuillez saisir un compte.'; the second is for the password with a lock icon and the placeholder text 'Veuillez saisir votre mot de passe.'. A blue 'SE CONNECTER' button is positioned below these fields.

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